

# Reichert 7CR

Auto Tonometer + Corneal Response Technology™

The Glaucoma Tonometer.



# **Reichert 7CR** Auto Tonometer + Corneal Response Technology™

IOP is the only modifiable risk factor for glaucoma and, as such, accurate tonometry is of critical importance in the diagnosis and management of this disease. It is widely recognized that the current gold standard for measuring pressure, the Goldmann Applanation Tonometer (GAT), has considerable flaws. Goldmann, and other current tonometers, are significantly influenced by corneal properties such as rigidity and thickness. Though it has been shown that GAT provides clinically reliable measurements in eyes with "average" corneas, we now know that many corneas vary more significantly from average than previously thought.

Numerous Studies have shown that thicker corneas overstate IOP values and thinner corneas understate IOP. However, thickness does not adequately describe corneal biomechanical properties, making the use of CCT-based correction algorithms questionable at best. In fact, researchers have concluded that correcting IOP based on CCT can lead not only to errors in the magnitude, but also in the direction of the adjustment.

"Correction nomograms that adjust GAT IOP based solely on CCT are neither valid nor useful in individual patients."


- Pg 18. Robert N. Weinreb, James D. Brandt, David Garway-Heath and Felipe Medeiros. World Glaucoma Association on Intraocular Pressure; Consensus Series 4; May 5, 2007

## **Reichert's 7CR** Auto Tonometer + Corneal Response Technology

takes corneal biomechanical properties into consideration, providing Corneal Compensated IOP (IOPcc) - a pressure measurement that is significantly less affected by the cornea than other methods of tonometry. The instrument's unique ability to simultaneously provide a Goldmanncorrelated IOP measurement (IOPq) and IOPcc, provides clinicians with a better understanding of patient tonometry values, enhancing clinicians ability to make critical diagnosis and treatment decisions.

#### Reichert 7CR Features

- IOPg (Goldmann Correlated IOP) and IOPcc (Corneal Compensated IOP) clearly displayed on color LCD
  - Best Signal Value (BSV) is automatically selected
- Easy-to-use, touch-screen user interface requires minimal operator training
- Exclusive Reichert technology ensures the softest possible air puff
  - Real-time applanation detection system provides a custom air puff for every measurement
- Simple and fast patient positioning No chinrest, joystick, or elevation controls
  - Patients simply lean against convenient left/right sliding forehead rest
  - Obvious fixation cues for fast alignment
  - Alignment and measurement is completely automated
- Triple-measurement mode delivers three consecutive IOP measurements with a single button press
- Internal printer and Electronic data transfer via USB port
- Reichert's Corneal Response Technology features patent US 7, 481, 767 B2

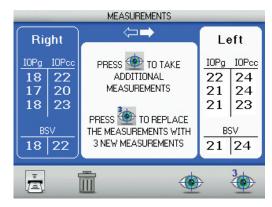

Reichert 7CR simultaneously displays IOPg and IOPcc

### Reichert 7CR is more accurate in:

- Normal Tension Glaucoma patients
- Primary Open Angle Glaucoma patients
- Post-LASIK and refractive surgery patients
- Fuchs' / edematous patients
- Keratoconus patients
- Patients with thick, thin, or biomechanically atypical corneas

# **Specifications**

Catalog Number	16060
Dimensions	19.75" h x 10.5" w x 14" d (50.2 x 26.7 x 35.6 cm)
Weight	23 lbs (10.4 kg)
Voltage	100/240 VAC
Frequency	50/60 Hz
Measurement Range	7 - 60 mmHg (ISO 8612 tonometer standard)





Reichert, Inc.

3362 Walden Avenue, Depew, NY 14043 Toll Free: 888-849-8955 • Tel: (716) 686-4500 Fax: (716) 686-4545 • www.reichert.com

ISO 9001 Certified

Reichert products are designed and manufactured under quality processes meeting ISO 9001 requirements.

Proudly Made in the USA.

#### **Call your authorized Reichert Ophthalmic Instruments Distributor today.**

**Important Notice**For your own protection, only purchase your Reichert Ophthalmic Instruments from and have them serviced by an Authorized Reichert Distributor.